★ ★ ★ ●
TRISTAR

WG - 2417

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Technical specification:

- | | |
|--|---------------------------------------|
| 1. Maxi measuring capacity: 180kg/396lb/19st | 2. Mini measuring capacity: 3kg |
| 3. Weight tolerance: +/-1% | 4. Weight graduation: 0.1kg/0.2 lb |
| 5. Weight unit: kg/lb /st | 6. Fat measuring range: 5%~65% |
| 7. Fat percentage graduation: 0.1% | 8. Water measuring range: 20%~72% |
| 9. Water percentage graduation: 0.1% | 10. Muscle measuring range: 20%~70% |
| 11. Muscle percentage graduation: 0.1% | 12. Calorie graduation: 1 kcal |
| 13. Memory units: 1-12 | 14. Personal height: 100~220cm |
| 15. Age range: 10~99 years | 16. Overloaded indication: 0-Ld |
| 17. Fat / water / muscle error indication: "Err" | 18. No saved data indication: "-----" |
| 19. Low battery indication: "Lo" | |

Directions for use


- Measuring weight

(1)

Put the scale on a level floor and stand still on it. The scale will turn on automatically and display the result.

You can also turn on the scale by pressing the button  and the weight result will flash two times before the scale turns off automatically (The scale will turn off automatically without operation within about 30 seconds when "0.0" appears).

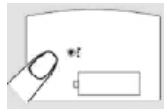(2)

St
Lb
Kg


Press the button at the back panel to change weight units when "0.0" shows on the display.

- Setting personal data


(1)


Press the button  to turn on the scale. Press and hold the button  for three seconds when "0.0" appears on the display, then personal memory unit will start flashing on the display, whose data can be adjusted.

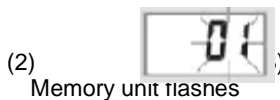
Memory unit flashes

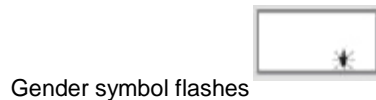
Gender symbol flashes

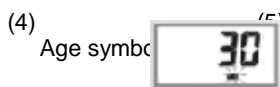Age symbol

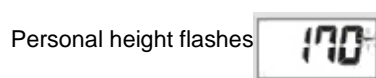Personal height flashes

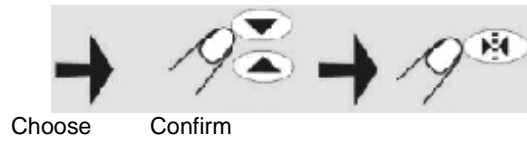(6) The memory unit appears on the display when the personal data are entered fully. Repeat the steps (2)~(5) to modify data or enter data for another user. The scale will shut off automatically without operation when the memory unit flashes 5 times and the incomplete setting will not be saved (The initial/default data are: memory unit "1", weight unit (kg), male, age 30, personal height 170cm).

Measuring body weight and analyzing fat, water and muscle percentage

Press the button " " to turn on the scale. Press the button " " to choose memory unit when "0.0" appears on the display, then press the button " " to enter. The personal data saved in the chosen memory unit will appear on the display and "0.0" appears subsequently.

(2)

Stand still on the scale when "0.0" appears, the weight value will be locked and the symbol appears circularly to indicate that the scale is measuring body fat, water and muscle, and is providing Calorie value for the person who should intake per day. The result will appear on the display alternately for 2 times before the scale turns off automatically.

(3) When measuring body fat, water and muscle, make sure that the feet are bare and clean and fully in touch with the pole on the surface to avoid composition percentage error. The scale will shut off automatically after the symbol "Err" flashes.

Precautions

Keep feet bare and clean

Do not pile up sundries

Do not immerse in water

Do not use or store

Do not whack

Avoid overloading

in damp environment

Standard Body Fat Percentage Distribution Table (%)

Age	Female %				Male %			
	Thin	Normal	High	Over	Thin	Normal	High	Over
15 – 39	14 – 20	21 – 25	26 – 31	32 – 38	11 – 17	18 – 23	24 – 29	30 – 36
40 – 55	15 – 21	22 – 26	27 – 32	33 – 38	12 – 18	19 – 24	25 – 30	31 – 37
56 – 70	16 – 22	23 – 27	28 – 33	34 – 38	13 – 19	20 – 25	26 – 31	32 – 38

Standard Body Water Percentage Distribution Table (%)

Age	Female %			Male %		
	Thin	Normal	High	Thin	Normal	High
10 – 15	< 57	57 – 67	> 67	< 58	58 – 72	> 72
16 – 30	< 47	47 – 57	> 57	< 53	53 – 67	> 67
31 – 60	< 42	42 – 52	> 52	< 47	47 – 61	> 61
61 – 80	< 37	37 – 47	> 47	< 42	42 – 56	> 56

Standard Body Muscle Percentage (For reference only)

Age	Female %	Male %
15 - 70	> 34	> 40%

WARNING: people with pacemakers or other electronic body implants should not use this product since the body resistance measuring method may interfere with its proper operation.

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Guidelines for protection of the environment

This appliance should not be put into the domestic garbage at the end of its useful life, but must be disposed of at a central point for recycling of electric and electronic domestic appliances. This symbol on appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling used domestic appliances you contribute an important push to the protection of our environment. Ask your local authorities for information regarding the point of recollection..

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